

How to don gloves

1



Take out a glove from its original box

2



Touch only a restricted surface of the glove corresponding to the wrist (at the top edge of the cuff)

3



Don the first glove

4



Take the second glove with the bare hand and touch only the restricted surface of glove corresponding to the wrist

5



To avoid touching the skin of the forearm with the gloved hand, turn the external surface of the glove to be donned on the folded fingers of the gloved hand, thus permitting to glove the second hand

6



Once gloved, hands should not touch anything that is not defined by indications and conditions for glove use

How to remove gloves

1




Pinch one glove near the wrist to remove it, without touching the skin of the forearm, and peel away from the hand, thus allowing the glove to turn inside out

2



Hold the removed glove in the gloved hand and slide the fingers of the ungloved hand inside between the glove and the wrist. Remove second glove by rolling it down the hand and fold into the first glove

3

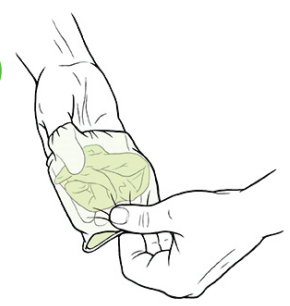

Discard the removed gloves

4

Then, perform hand hygiene by rubbing with an alcohol-based handrub or by washing with soap and water